


Accessories

	Master platform set Art.-No. 72203911 Supplied: individually		No-Plast Art.-No. 72203912 Supplied: 1 pack of 2
	Calibration key Art.-No. 72203913 Supplied: calibration key with adjustment plate		Face bow support Art.-No. 72203930 Supplied: individually
	Lower adjustment plate for calibration Art.-No. 72203914 Supplied: individually		Denar facebow adapter Art.-No. 7220 3940 Supplied: individually
	Instrument holder with Flat setting up template Art.-No. 72203915 Supplied: complete		Face-bow adapter for: • SAM Axioquick AX • Kavo Protar Arcus • Girrbaach Artex® face-bow • Ivoclar UTS transfer-bow Art.-No. 72203921 Supplied: 1 pack of 2 (L/R)
	Flat setting up template Art.-No. 72203920 Supplied: individually		

Spare Parts

Bi-conical axle 13° - 17° Art.-No. 72203900 Supplied: 1 pack of 2 (R/L)	Incisal guide Art.-No. 72203901 anodised AL 0° 72203902 anodised AL 15° standard 72203903 anodised AL 30° Supplied: individually	Incisal pointer Art.-No. 72203904 Supplied: individually
Plastic mounting plates Art.-No. 72203905 Supplied: 1 pack of 10	Metal mounting plates Art.-No. 72203906 Supplied: 1 pack of 2	Retention discs Art.-No. 72203907 Supplied: 1 pack of 100
Plastering covers Art.-No. 72203908 Supplied: 1 pack of 2	Labels Art.-No. 72203909 Supplied: 1 pack of 2	O ring seals Art.-No. 72203910 Supplied: 1 Pack of 2
Vertical guide rod Art.-No. 72203922 Supplied: individually	Vertical guide rod holder Art.-No. 72203923 Supplied: individually	Magnet Art.-No. 72203924 Supplied: 1 Pack of 2
Conical nut Art.-No. 72203925 Supplied: 2 pieces	Knurled thumb screw A Art.-No. 72203926 Supplied: 2 pieces	Knurled thumb screw B Art.-No. 72203927 Supplied: 2 pieces
Centric quick-lock incl. screw Art.-No. 72203928 Supplied: individually	Supporting rod opening stop Art.-No. 72203929 Supplied: individually	Magnetic retention plate Art.-No. 72300510 Supplied: 2 pieces

Candolor Articulator CA II



The Prosthetic Articulator

Candolor Articulator CA II

The Prosthetic Articulator

Articulators play a decisive role in the success of our restorations because occlusion and articulation always have to be considered, even for the smallest dental prosthesis.



Stable

A sturdy construction in light-weight aluminium, the striking features of the Candolor Articulator are its clean lines, compact size and simple handling.

This prosthetic articulator is a semi-adjustable articulator designed in accordance with Camper's plane.



Simple

The articulator is equipped with a split-cast system with strong magnets for fast, easy fixing of the mounting plates. An important design feature is the calibration key for achieving commonality between several Candolor articulators, which enables models to be interchanged at any time.

The split-cast recesses incorporated in the upper and lower articulator plates enable the model to be directly fixed with plaster. This makes the base suitable for grinding to give the finished work a neater appearance.

Universal

The Candolor prosthetic articulator features a transfer system which makes it compatible with joint and skull-related face-bows (see list of accessories).

Functional


The mechanical joint takes the form of a guide edge which supports a bi-conical axle.

The bi-conical design emulates the shape of the condyle to prevent straightline movements in the case of lateral and transverse displacements, thus ensuring adherence to physiological principles.

The joint axle is opened and closed easily by means of a precise centric quick-lock and retrusion lock.

The horizontal condyle path inclination can be adjusted from 0° – 60° by rotating the wheel.

Features / Benefits

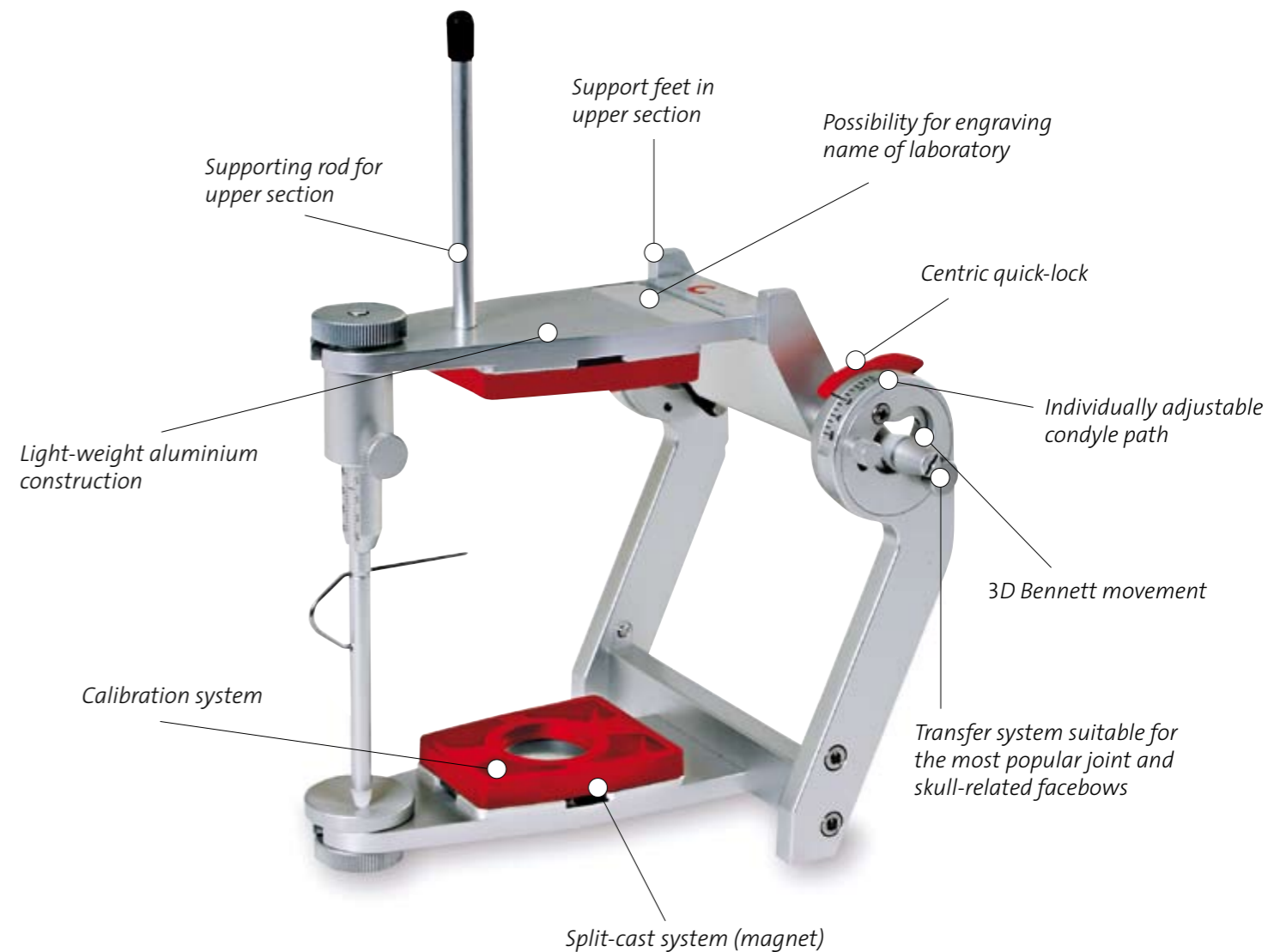

2-year warranty against material failure

Supplied complete

- 1 x Candolor Articulator CA II
- 1 x Incisal guide 15°
- 1 x Incisal pointer
- 2 x Retention discs
- 2 x Mounting plates
- 2 x Labels
- 1 x Supporting rod for upper section (opening stop)
- 2 x Plastering covers
- 2 x Spare O ring seals
- 2 x Magnetic retention plates

Specifications

Height: 150 mm
 Inside height: 105 mm (with split-cast)
 Width: 145 mm
 Depth: 160 mm
 Weight: 690 g

- Bonwill's triangle 110 mm
- Balkwill angle 20°
- Retrusive path 1.5 mm
- Immediate side shift stop: fixed 2.5 mm
- Condyle path inclination adjustable 0° – 60°
- Material: anodised aluminium
- Colour: silver